


Logitech Mouse G700 Setup


Setup of Logitech Wireless Gaming Mouse G700 for Digital Microscope

Connect and setup the mouse

- **Unpack and check the package contents**

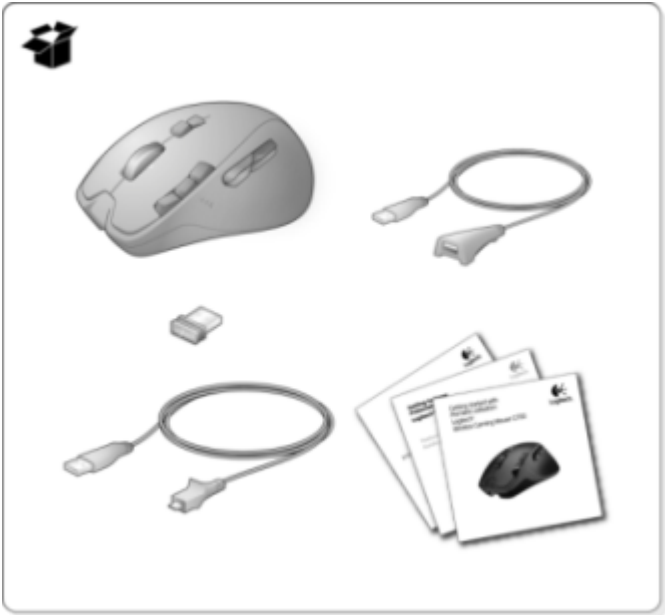


- **Activate G700**

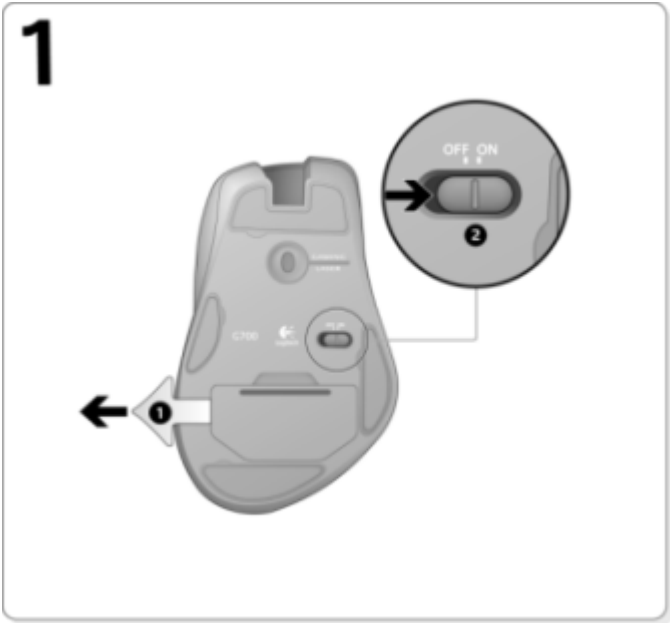

- **Connect the receiver. .**

2

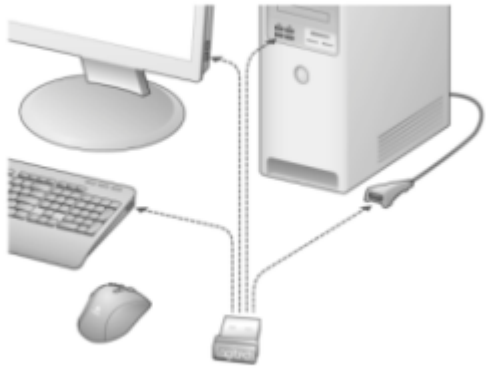

- Install Logitech Gaming Software from <http://www.logitech.com/en-us/support>


3

www.logitech.com/downloads



- Charge the mouse.

4

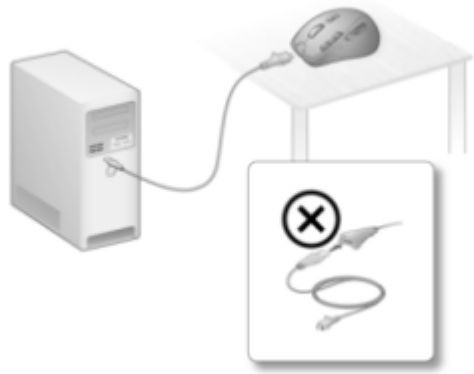

Create profile for Digital Microscope

- Start **Logitech Gaming Software**.



- Select **Automatic Game Detection** in the **Home** tab.




- Select **Customize Buttons** tab.



- Click **Import/Export** button (folder icon) and select **Import...** in the menu.



- Find **DM-Logitech-G700-Custom-Profile.xml** in DM installation folder and press "Open".

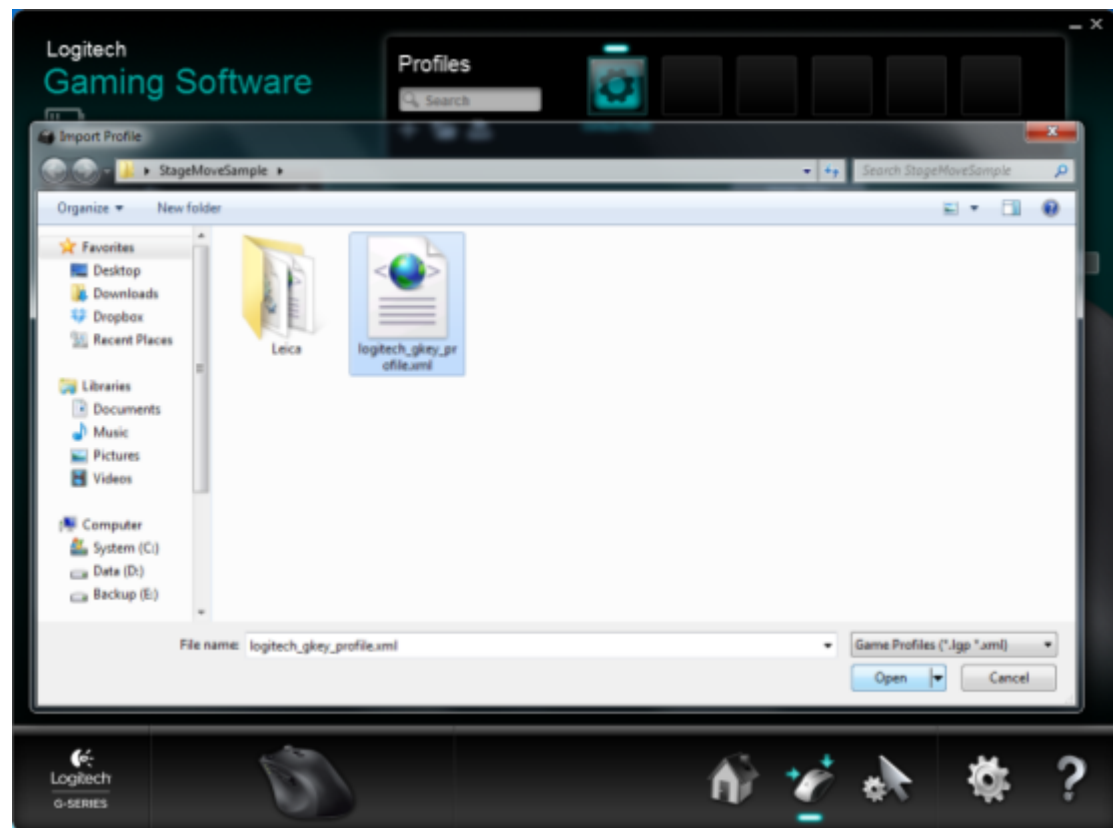


- In **Profiles** list there should be new **DM** profile, click it with right mouse button and select **Properties**.



- In Edit Profile window click **Select Game Executable**.



- Find DMClient.exe in DM application folder and click "Open".

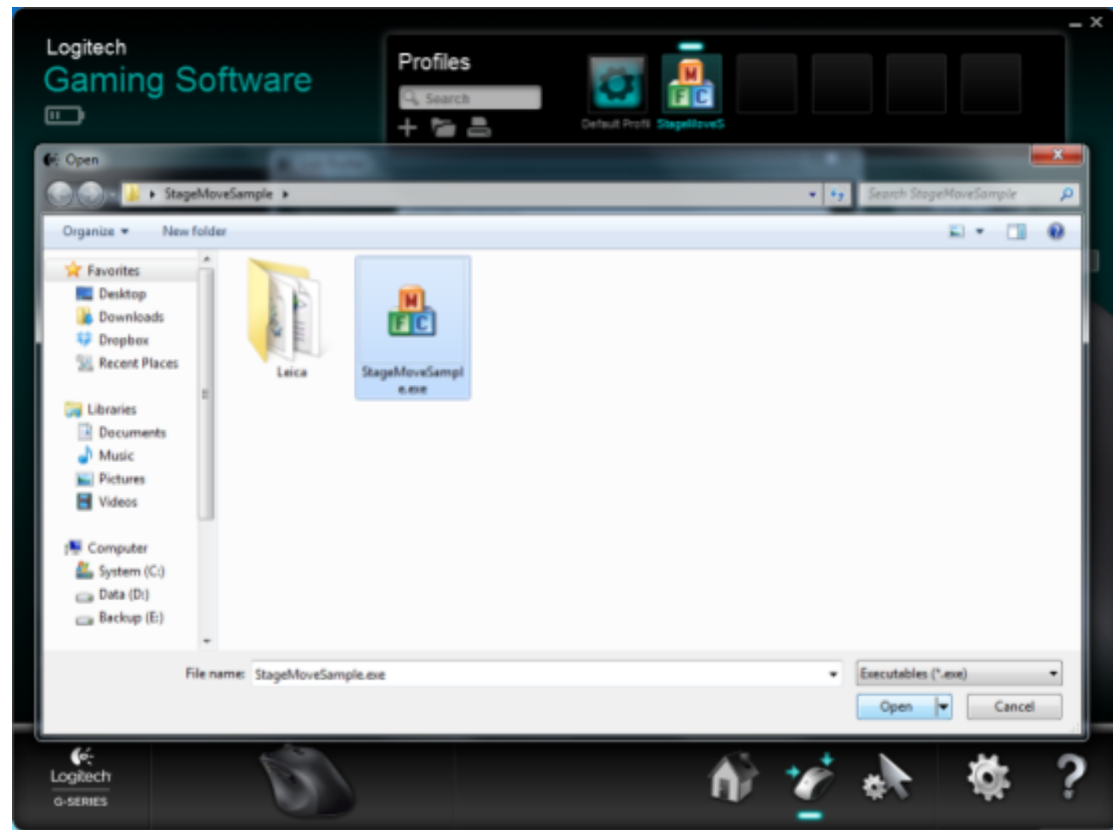

Profile installed is finished. Now it will be activated on each start of DM application.