

# DiBox 2.0 data for inclusions plotting


Darkfield and Brightfield images with recognized basis along with dmc model uploaded on Cutwise are required to use OSV Plotting.

To receive photo with recognized basis:

1. Switch to basis recognition mode in DiBox software Edit panel:




2. Load dmc model to dibox

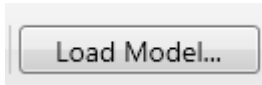

3. Recognize basis on ASET white (or ASET Black) photo



4. Copy to darkfield or brightfield photo

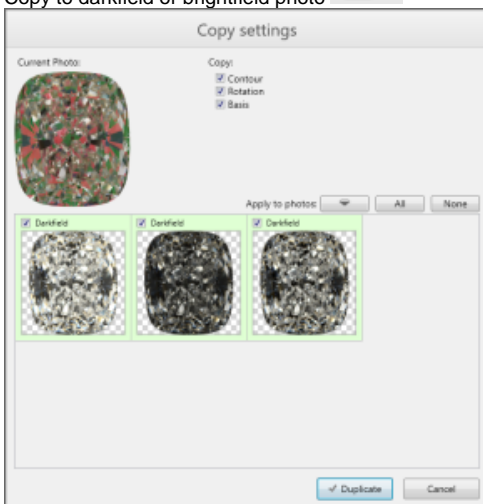

- 5.

Upload images with recognized basis along with .dmc file to Cutwise.

Start OSV Plotting and plot inclusions.